Instituto Tecnológico y de Estudios Superiores de Monterrey



BI2009B.300 Procesamiento de imágenes médicas para el diagnóstico

Profesor: José Gerardo Tamez Peña Profesora: Martha Rebeca Canales Fiscal

Entregable 1

MRI Activity

Marcela Enriquez López | A01570502

Natalia Verónica Flores Del Río | A01570472

Ana Lucía Soria Cardona | A00827565

Mei Li Luisa Cham Pérez | A01139386

Graciela Alejandra Rincón López | A00827270

MRI activity.

Open the simulator in https://phet.colorado.edu/en/simulations/mri

1. NMR

Place yourself in the tab *Simplified NMR*. Remember that the Larmour relationship relates the magnetic field to the resonant frequency:

$$v_L = \frac{\gamma}{2\pi} B_0$$

where v_L is the resonant frequency, $\frac{\gamma}{2\pi}$ is the gyromagnetic ratio and B_0 is the magnetic field.

Table 1. Gyromagnetic ratios for different nuclei.

Nuclei	Gyromagnetic ratio	Nuclei	Gyromagnetic ratio
1H	42,58	65Cu	12,09
7Li	16,55	75As	7,291
9Be	5,984	77Se	8,118
11B	13,66	81Br	11,50
13C	10,71	87Rb	13,93
15N	4,314	93Nb	10,41
170	5,772	117Sn	15,17
19F	40.05	121Sb	10,19
23Na	11,42	1271	8,518
27AI	11,09	133Cs	5,584
29Si	8,458	195Pt	9,153
31P	17,24	199Hg	7,590
35CI	4,172	203TI	24,33
51V	11,19	207Pb	8,907
55Mn	10,50	209Bi	6,841
59Co	10,05		


Use the Larmour relationship and the gyromagnetic ratios of various nuclei shown in table 1 to complete table 2. Check your results in the simulation by setting the appropriate frequencies and main magnet field, take a screenshot of the nuclei emitting energy to


include in the report. Try to find the last nuclei (????) by playing with the simulation and register the frequency at two different magnetic fields.

Table 2. Different settings to achieve energy emission.

Nuclei	Magnetic Field	Resonant Frequency	Magnetic Field	Resonant Frequency
Hydrogen	0.75	31.935	2.5	106.45
Nitrogen	2.5	10.785	-	-
Sodium	1.515	17.3	2.75	31.405
Carbon-13	1.75	18.742	2.5	26.775
Oxygen	1.999	11.54	3.0	17.316
Sulfur	3.0	25.375	-	-
????	1.75	22.0	3	36.4

The last nuclei is Copper

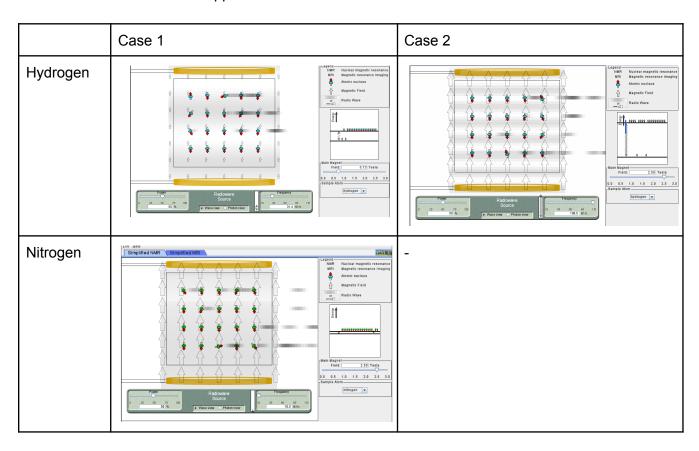

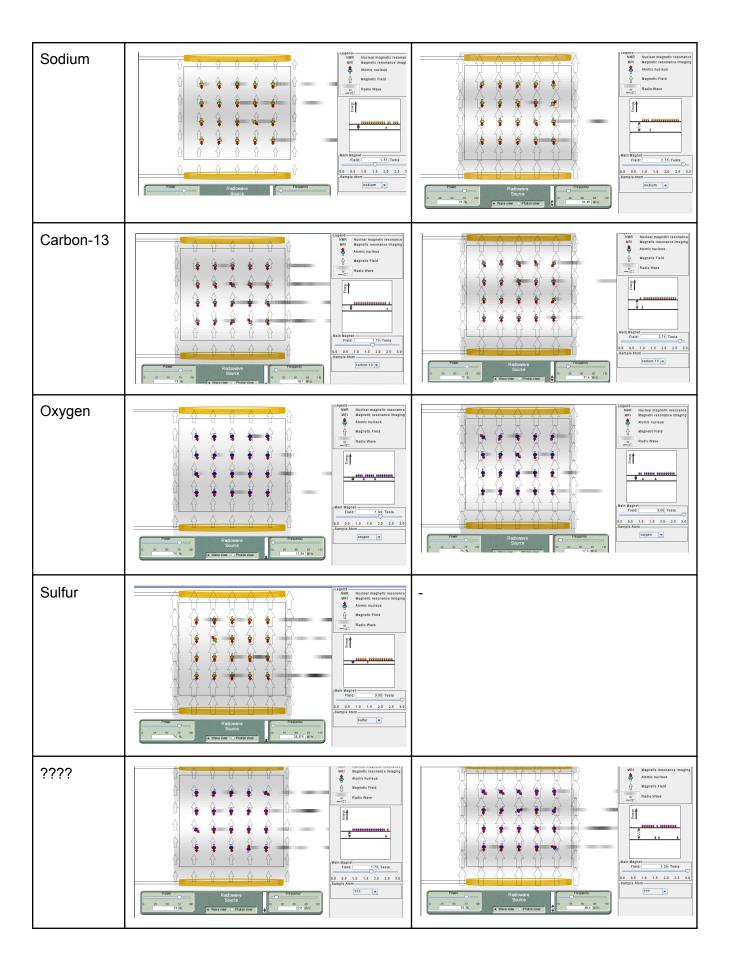


2. MRI

Move to the Simplified MRI tab


a. Set the main magnet field to 1.0 Tesla, leave the gradient magnets in zeros, activate only show head, and show magnetic field (be sure that show atomic nuclei is deactivated), set the frequency in 43 MHz. Finally set the power to 50% and observe the flow and distribution of the emissions. After a while observing the emissions, click on add tumor, wait for around 7 seconds so the flow distribution stabilizes, look at how the emission changed and try to guess where the tumor is located.

Explain how the emission allowed you to find the correct location:

Since the tumor causes an obstruction in the electromagnetic flow its repercussions on the flow can be key to determine the presence of a tumor, due to the fact that when the flow hits the tumor its flow is disturbed. This can be seen throughout the first moments that the flow is severy disrupted, specially in the areas surrounding the tumor. Once it stabilizes the flow surrounding the tumor is still chaotic, yet the flow on the other sections of the mri is more stable. The detection can be found due to the disruptions on the flow.

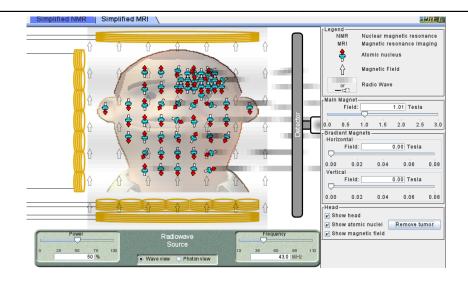

b. Play with the main magnet field, frequency, and gradient magnets (both, horizontal and vertical) to try to obtain an emission focused mainly in the zone of the tumor (register your best guess, it doesn't need to be perfect). Answer the following questions.

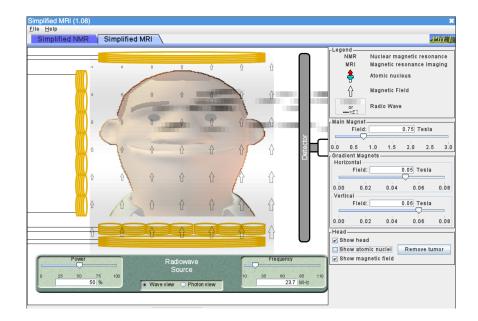

Best guess: main magnet: <u>0.75 Tesla</u> horizontal gradient: <u>0.05 Tesla</u> vertical gradient: <u>0.06 Tesla</u> frequency: <u>23.7 MHz</u>

What happens when the horizontal gradient increases its magnitude? How does it affect the emissions? When the horizontal gradient increases the concentration of the emissions in the area decreases

What about vertical gradient? When the vertical gradient increases the emission increases in the area of the tumor and the emissions are more constant.